GRUNDIG

Sonic Toothbrush

User Manual

TB 8730

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01M-GMS2300-0321-06

DEUTSCH 06-21

ENGLISH 22-32


TÜRKÇE 33-49











ESPAÑOL 50-62

FRANÇAIS 63-75

HRVATSKI 76-88

POLSKI 89-101

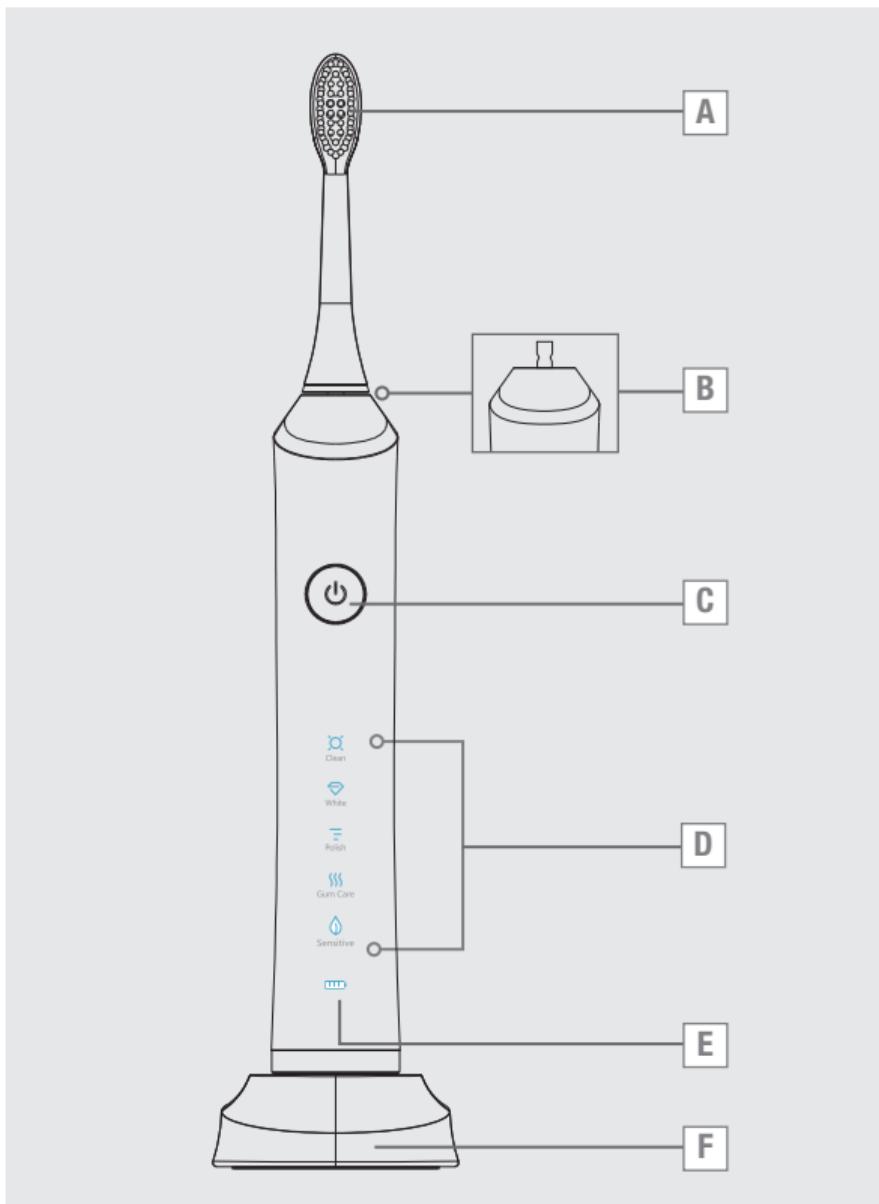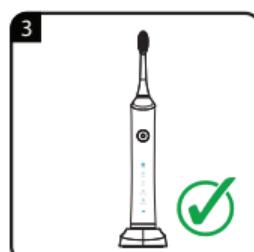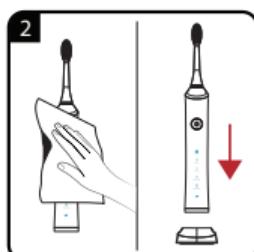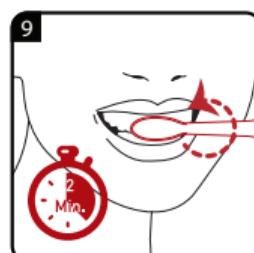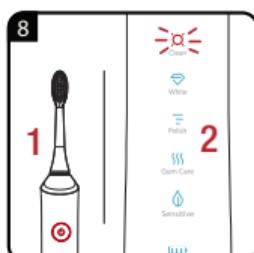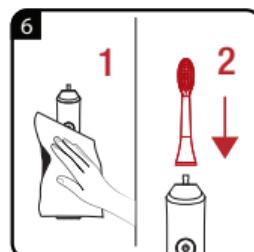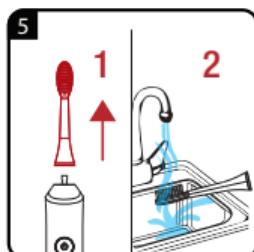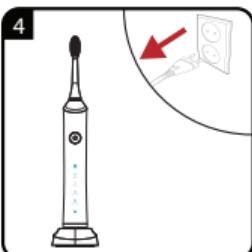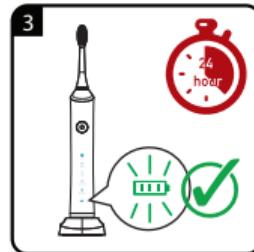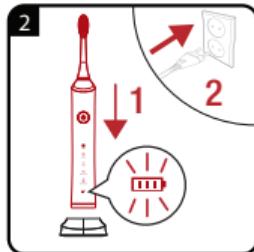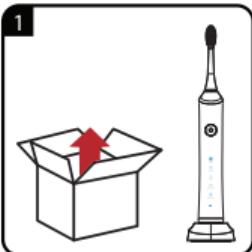		
DE	BETRIEB	REINIGUNG UND PFLEGE
EN	OPERATION	CLEANING AND CARE
TR	KULLANIM	TEMİZLEME VE BAKIM
ES	FUNCIONAMIENTO	LIMPIEZA Y CUIDADOS
FR	FONCTIONNEMENT	NETTOYAGE ET ENTRETIEN
HR	RAD	ČIŠĆENJE I ODRŽAVANJE
PL	OBSŁUGA	CZYSZCZENIE I KONSERWACJA

SICHERHEIT

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Garantiebedingungen - Grundig - 01/2021

SAFETY

Please read this instruction manual thoroughly prior to using this appliance! Follow all safety instructions in order to avoid damage due to improper use!

This appliance is designed solely for cleaning and caring for human teeth. All other uses are expressly prohibited.

- Keep the instruction manual for future use. Should a third party be given the appliance, please ensure the instruction manual is included.
- Do not use if the appliance, attachments or the power cord on the charging unit is visibly damaged.
- If the power cord from the charging unit is damaged, it is not possible to repair it. Please dispose the set according to environmental rules.

SAFETY

- When the appliance is being charged, on no account should it or the charging unit come into contact with water or any other liquid.
- Keep the appliance out of direct sunlight, excessive heat and open flames.
- Do not drop the appliance.
- Keep the appliance out of the reach of children.
- Toothbrushes can be used by children, by people who have reduced physical, sensory or mental capabilities and by people who lack experience or knowledge on the condition that they are supervised or instructed about the safe use of the product and if they understand the possible dangers involved.

SAFETY

Cleaning and user maintenance should not be performed by children without supervision. Children should not play with the appliance.

- Do not open the appliance under any circumstances. No warranty claims are accepted for damage caused by improper handling.
- This product is intended to be used with water or toothpaste only. Do not use it with anything else.
- Use only the parts and adaptor provided with the product.
- Make sure that you unplug the product before cleaning and maintenance.

Dear Customer,
congratulations on purchasing the GRUNDIG Sonic Toothbrush TB 8730.

Please read the following user guide carefully to ensure that you have many years of enjoyment from your quality Grundig product at home.

Responsible trading!

Both internally and for suppliers, GRUNDIG insists on contractually guaranteed social working conditions with fair wages, efficient use of raw materials with continuous reduction in waste of several tonnes of plastic each year - and at least 5 years availability for all accessories.

For a future worth living.

Grundig.

Controls and parts

See the figure on page 3.

- [A] Brush head attachment
- [B] Brush head holder
- [C] On/Off and mode switch
- [D] Integrated display for different settings
- [E] Indicator LED for charging status
- [F] Charging unit

Using in the bathroom

Your toothbrush has a waterproof handle, is electrically safe and can be used in the bathroom without hesitation.

Oral hygiene

Careful oral hygiene is essential for maintaining healthy, beautiful teeth.

Daily intensive oral hygiene can prevent tooth decay, gum disease, inflammation of the gums and bad breath.

Brushing time

The recommended time for brushing is 2 minutes. The timer then stops the function automatically.

Settings

With 5 different settings:

- **Clean:** 2-minutes cleaning mode
- **White:** 2-minutes cleaning mode with additional 30-seconds brightening and polishing of front teeth
- **Polish:** 1-minute fast mode to brighten and polish front teeth
- **Gum Care:** 2-minutes cleaning mode with additional 1-minute gum-and teeth root massage, 15-seconds timer reminds changing direction of brush
- **Sensitive:** 2-minutes extra sensible cleaning mode for protection of teeth and gum.

Charging

Check that the mains voltage on the type plate of the charging unit (on the bottom) corresponds to your local mains supply.

Notes

- Charge the batteries fully before you use the appliance for the first time. It takes approximately 24 hours (12 hours charging time afterwards) to charge them initially. When the appliance is fully charged, the charge LED  lights up green.
- The ambient temperature should be between 5 °C and 40 °C.

- When the batteries are fully charged, they provide approximately 75 days of operation. The charging process ends automatically once the batteries are fully charged. We recommend keeping the appliance in the charging unit so that it is available fully charged at all times.

DENTAL CARE

Notes

- Rinse your mouth out with water before brushing your teeth.
- The chosen programm name lights on the display.
- Integrated timer indicates optimum usage time of 2 minutes in 30 seconds intervals.
- The toothbrush works best if you apply slight pressure when brushing and continue with the circular movements at the same time.
- In some cases you might experience a tickling or tingling sensation when using an electric sonic toothbrush for the first time. This sensation gradually subsides once you have got used to the sonic toothbrush.

Notes

- We also recommend using a mouth wash.
- Consult your dentist if you feel any sensitivity while using the appliance for the first time.

Cleaning and care

Cleaning

The sonic electric toothbrush removes plaque more precisely and effectively than conventional manual toothbrushes due to the high speed of the brushing movements (up to 48.000 rpm) on the surface of the teeth.

Discolouration of the teeth from coffee, tea or tobacco can be removed in a natural way.

Teeth, tongue and gums can be cleaned and cared for carefully.

Warning

- Never use gasoline, solvent, abrasive cleaning agents, metal objects or hard brushes to clean the appliance.

Note

- Replace the brush head when the bristles become bent and distorted and the wear indicator becomes brighter at the side of the bristles. We recommend you to do this every 3 months.
- Rinse the brushhead after each usage.

Storage

- If you do not plan to use the appliance for a long period of time, please store it carefully.
- Make sure that it is unplugged and is completely dry.
- Store the appliance in a cool, dry place.
- Make sure the appliance is kept out of the reach of children.

Compliance with the WEEE Directive and Disposing of the Waste Product:

This product complies with EU WEEE Directive (2012/19/EU).

This product bears a classification symbol for waste electrical and electronic equipment (WEEE).

This symbol indicates that this product shall not be disposed with other household wastes at the end of its service life. Used device must be returned to official collection point for recycling of electrical and electronic devices. To find these collection systems please contact to your local authorities or retailer where the product was purchased. Each household performs important role in recovering and recycling of old appliance.

INFORMATION

Appropriate disposal of used appliance helps prevent potential negative consequences for the environment and human health.

Information on Waste Batteries

This symbol on rechargeable batteries/batteries or on the packaging indicates that the rechargeable battery/battery may not be disposed of with regular household waste.

Always dispose of used rechargeable batteries/batteries in accordance with local environmental regulations. Make enquiries about the applicable disposal regulations where you live.

Compliance with RoHS Directive

The product you have purchased complies with EU RoHS Directive (2011/65/EU). It does not contain harmful and prohibited materials specified in the Directive.

Package information

Packaging materials of the product are manufactured from recyclable materials in accordance with our National Environment Regulations. Do not dispose of the packaging materials together with the domestic or other wastes.

Take them to the packaging material collection points designated by the local authorities.

Electrical insulation class: II

Technical data

Power supply

Charging unit: 100-240 V~,
50-60 Hz

Output rating: 4V  3W

Batteries: Li-ion (Lithium-ion rechargeable battery)

Protection class: (IPX7)

Design and specifications subject to change without notice.

Please see the cover page of the user manual for initial use, general use and cleaning sections.

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Beko Grundig Deutschland GmbH